


This two-part Magnetic System conveniently reduces the cycle time to process microstructure materials specimens. The two part magnetic based system consists of a MAGNETIC PSA (pressure sensitive adhesive) backed DISC, the PSA is directly attachable to existing platens. The MAGNETIC DISC receives a (second) METALLIC DISC, the versatile, non-corrosive METALLIC DISC accepts a complete variety of PSA abrasive discs, polishing cloths and pads, and all available diamond PSA products used during the plan grinding, fine grinding and polishing process.

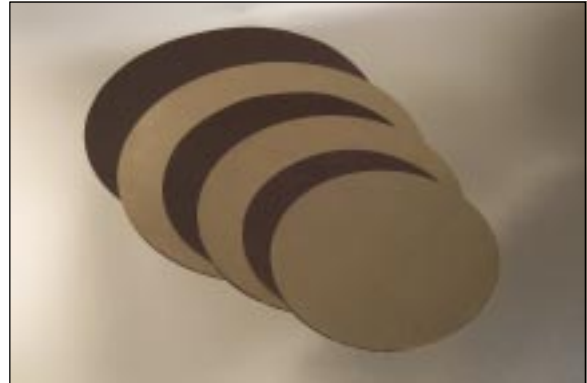

CATALOG NUMBER	DESCRIPTION	UNIT
EXTEC Magnetic Disc System		
16008	Magnetic Disc 8" (203 mm) PSA	Each
16010	Magnetic Disc 250 mm (9.85") PSA	Each
16012	Magnetic Disc 12" (305 mm) PSA	Each
16009	Metallic Disc 8" (203 mm)	Each
16011	Metallic Disc 250 mm (9.85")	Each
16013	Metallic Disc 12" (305 mm)	Each